

Takara Bio USA

# Make Your Own “Mate & Plate™” Library System User Manual

Cat. No. 630490  
(122321)

---

**Takara Bio USA, Inc.**

2560 Orchard Parkway, San Jose, CA 95131, USA

U.S. Technical Support: [technical\\_support@takarabio.com](mailto:technical_support@takarabio.com)

United States/Canada  
800.662.2566

Asia Pacific  
+1.650.919.7300

Europe  
+33.(0)1.3904.6880

Japan  
+81.(0)77.565.6999

Page 1 of 25

## Table of Contents

I.	Introduction & Protocol Overview.....	4
	A. About this Manual.....	4
	B. Matchmaker Library Construction Overview .....	5
II.	List of Components .....	5
III.	Additional Materials Required.....	6
IV.	List of Abbreviations.....	7
V.	Control Experiments .....	7
	A. General Considerations .....	7
	B. Protocol: Control for Homologous Recombination-Mediated Cloning Protocol .....	7
VI.	Generating the cDNA for Your Library .....	9
	A. Protocol: First Strand cDNA Synthesis.....	9
	B. Protocol: Amplify cDNA Using Long Distance PCR (LD-PCR) .....	11
	C. Protocol: Purify ds cDNA with CHROMA SPIN+TE-400 Columns .....	12
VII.	Construction of a Two-Hybrid Library .....	14
	A. General Considerations for Two-Hybrid Library Construction & Screening .....	14
	B. Protocol: Creating a Mate & Plate Library .....	14
VIII.	References .....	16
	Appendix A. Troubleshooting Guide.....	17
	Appendix B. Library Titering .....	19
	A. General Considerations .....	19
	B. Library Titering .....	20
	Appendix C. Yeast Growth Media & Supplements .....	20
	A. Ready-to-go Media Pouches Available from Takara Bio .....	20
	B. General Media Preparation Instructions.....	22
	C. Freezing Medium .....	22
	D. Kanamycin Supplement to Yeast Media .....	22
	Appendix D. Host Strain Information.....	22
	Appendix E. Plasmid Information .....	23
	Appendix F. SMART Technology Overview .....	24
	A. SMART Technology .....	24
	B. Mechanism of cDNA Synthesis .....	24

## Table of Figures

Figure 1. Use yeast biology to construct a Mate & Plate library .....	4
Figure 2. Mate & Plate Library construction overview .....	9
Figure 3. High-quality cDNA generated using SMART cDNA synthesis .....	12
Figure 4. CHROMA SPIN column and collection tubes .....	13
Figure 5. Restriction Map and Cloning Sites of pGADT7-Rec .....	23

## Table of Tables

Table 1. Make Your Own "Mate & Plate" Library System components .....	5
Table 2. Relationship between amount of RNA and optimal number of thermal cycles .....	11
Table 3. Troubleshooting guide .....	17
Table 4. Components of Yeast Media Set 2 & Yeast Media Set 2 Plus .....	21
Table 5. Individual Yeast Media Pouches for Matchmaker Gold Protocols.....	21
Table 6. Additional Media and Media Supplements Required for a Two-Hybrid Screen.....	22
Table 7. Yeast Host Strain Genotype.....	22
Table 8. Phenotype Testing on Various SD Media.....	22


## I. Introduction & Protocol Overview

Yeast two-hybrid systems are primarily used for screening a complete library of proteins (prey) for interaction with a specific protein of interest (bait). Since traditional library manufacture is time-consuming and labor-intensive, Takara Bio has developed the Mate & Plate Libraries, a series of ready-to-go libraries that require simple co-culturing of a library strain with your bait strain followed by selection on appropriate minimal medium.

Our **Make Your Own “Mate & Plate” Library System** (Cat. No. 630490) enables you to “make your own Mate & Plate Library” just the way we do it using a simple and highly efficient protocol.

Library construction takes place directly in our library yeast strain Y187, utilizing the highly potent homologous recombination machinery of *Saccharomyces cerevisiae* (Figure 1). There is no need for labor-intensive library cloning, amplification, and harvesting in *E. coli*. This system uses SMART® cDNA synthesis technology, which allows you to construct cDNA libraries from any tissue source starting with as little as 100 ng of total RNA.

Refer to Appendix D for Y186 host strain information, Appendix E for detailed plasmid information, and Appendix F for an overview of SMART technology.

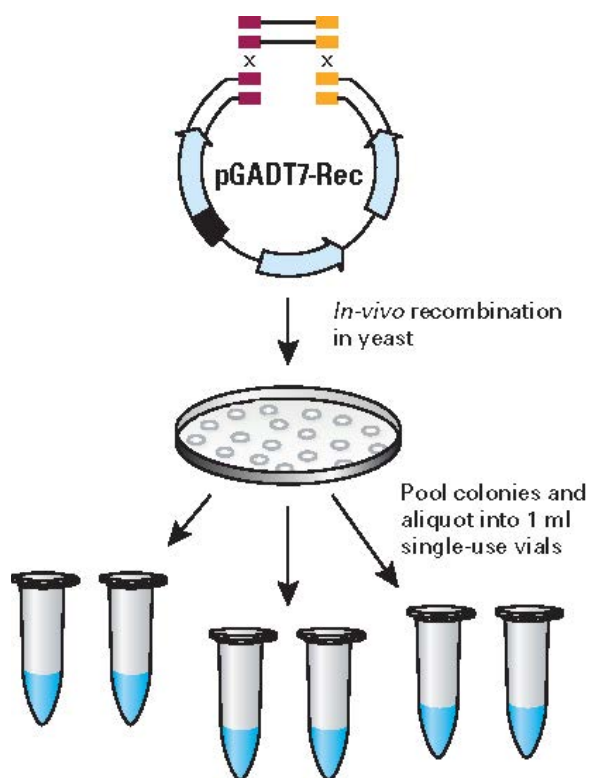

**Figure 1. Use yeast biology to construct a Mate & Plate library.** Mate & Plate libraries are made directly in yeast, via in vivo recombination between your cDNA and the Matchmaker prey vector pGADT7-Rec. First, use SMART cDNA synthesis technology to create a cDNA pool with end sequences homologous to the prey vector pGADT7-Rec. Then transform Y187 Yeast Strain and allow the yeast to perform a recombination step between the linear prey vector and the cDNA. Colonies are pooled, mixed, and aliquoted into multiple vials. Each single 1-ml vial can be used for a two-hybrid screen.

## A. About this Manual

This manual describes how to make your own library and divide it into aliquots in preparation for a yeast two-hybrid screen. For detailed procedures regarding the library screening process itself, please see our [Matchmaker® Gold Yeast Two-Hybrid System User Manual](http://www.takarabio.com) at [takarabio.com](http://www.takarabio.com). The Matchmaker Gold Yeast Two-Hybrid System is sold separately as Cat. No. 630489.

## B. Matchmaker Library Construction Overview

The library construction process consists of the following steps (see Figure 1):

1. Generate cDNA from RNA using SMART technology
2. Cotransform cDNA and linearized pGADT7-Rec vector into Yeast Strain Y187 and plate on SD/-Leu Agar medium
3. Wait 4–5 days, harvest yeast, and pool all colonies in rich broth medium.
4. Divide into aliquots and freeze.
5. Use a single 1-ml aliquot for each library screen (referenced in Section I.A).

## II. List of Components

Make Your Own “Mate & Plate” Library System (Cat. No. 630490) contains sufficient reagents to make five two-hybrid libraries.

**Table 1. Make Your Own “Mate & Plate” Library System components.**

<b>Make Your Own “Mate &amp; Plate” Library System</b>	<b>630490</b>
<b>Box 1 (Store at –20°C)</b>	
SMART MMLV RT (200 units/μl)*	10 μl
5X First-Strand Buffer	300 μl
DTT (100 mM)	165 μl
CDS III Primer (12 μM)	10 μl
CDS III/6 Primer (10 μM)	10 μl
RNase H (2 units/μl)	7 μl
5' PCR Primer (10 μM)	50 μl
3' PCR Primer (10 μM)	50 μl
Melting Solution	500 μl
SV40 Large T PCR Fragment (25 ng/μl)	20 μl
dNTP Mix (10 mM each dNTP)	50 μl
pGADT7-Rec AD Cloning Vector (Smal-linearized; 500 ng/μl)	25 μg
<b>Box 2 (Store at –20°C)</b>	
Yeastmaker™ Yeast Transformation System 2 (Box 1 of 2)†	
Yeastmaker Carrier DNA, denatured (10 mg/ml)	2 x 1 ml
pGBT9 (100 ng/μl; control plasmid)	20 μl
<b>Box 3 (Store at room temperature)</b>	
Yeastmaker™ Yeast Transformation System 2 (Box 2 of 2)†	
50% PEG	2 x 50 ml
1 M LiAc (10X)	50 ml
10X TE Buffer	50 ml
YPD Plus Liquid Medium	50 ml
<b>Box 4 (Store at –70°C)</b>	
SMART III Oligonucleotide (12 μM)	10 μl
Control Poly A <sup>+</sup> RNA (Mouse Liver; 1 μg/μl)	5 μl
<i>S. cerevisiae</i> Y187	0.5 ml
<b>Box 5 (Store at room temperature)</b>	
YPDA Broth (0.5 L)	1 pouch
YPDA with Agar (0.5 L)	1 pouch
SD/-Leu with Agar (0.5 L)	1 pouch
NaCl Solution (0.9%)	50 ml
Sodium Acetate (3 M; pH 4.8)	300 μl
Deionized H <sub>2</sub> O	500 μl
CHROMA SPIN™+TE-400 Columns	10 each

\*Also sold separately as Cat. No. 639522.

†Also sold separately as part of Cat. No. 630439.

### III. Additional Materials Required

- Advantage® 2 Polymerase Mix (100 rxns; Takara Bio, Cat. No. 639201)

#### Accessory Kits

- Yeastmaker Yeast Transformation System 2 (supplied with your system—see Section II; also sold separately as Takara Bio, Cat. No. 630439)
- Easy Yeast Plasmid Isolation Kit (50 preps; Takara Bio, Cat. No. 630467)

#### Tools for Plating Yeast

- Tools for plating yeast: includes a sterile glass rod and a bent Pasteur pipette or 5 mm glass beads for spreading cells on plates. (Use 5–7 beads per 100 mm plate, or 15–20 beads for a 150 cm plate)

#### Miscellaneous Reagents & Media

- 3 M sodium acetate, pH 5.3
- Ice-cold ethanol (95–100%)
- YPDA/25% Glycerol (Freezing Medium; see Appendix C, Section C)

#### Yeast Media Pouches Required for Library Construction (also see Appendix C)

- YPDA Broth (10 x 0.5 L; Takara Bio, Cat. No. 630306)
- YPDA with Agar (10 x 0.5 L; Takara Bio, Cat. No. 630307)
- SD/-Leu with Agar (10 x 0.5 L; Takara Bio, Cat. No. 630311)

#### Additional Yeast Media for Two-Hybrid Screening

(also see the [Matchmaker Gold Yeast Two-Hybrid System User Manual](#))

Table 4 (in Appendix C) lists the components of the **Yeast Media Set 2** (Takara Bio, Cat. No. 630494) and the **Yeast Media Set 2 Plus** (Takara Bio, Cat. No. 630495). These media sets contain a complete assortment of mixes for preparing eight specialized broth and agar media, designed for use with the Matchmaker Gold Yeast Two-Hybrid System, in convenient, “ready-mixed” foil pouches. The Yeast Media Set 2 Plus also contains the additional media supplements Aureobasidin A and X- $\alpha$ -Gal, which are required for a two-hybrid screen.

Table 5 (in Appendix C) contains information for purchasing each of the media mixes separately, in packs of 10 pouches, and Table 6 (in Appendix C) contains preparation instructions for all additional required media supplements and information for purchasing Aureobasidin A and X- $\alpha$ -Gal separately.

Additionally, the following should be considered when culturing yeast for a two-hybrid screen:

- **SD medium** (synthetically defined medium) is minimal media that is routinely used for culturing *S. cerevisiae*. SD base supplies everything that a yeast cell needs to survive (including carbon and nitrogen sources). Essential amino acids, which are added to SD base to create minimal medium, are already included premixed in Yeast Media Pouches from Takara Bio. The particular minimal medium that is chosen will determine which plasmids and/or activated reporters are selected for.
- **SD/-Leu dropout supplement** (either SD/-Leu with Agar (included in the kit; Takara Bio, Cat. No. 630311) or DO Supplement –Leu (sold separately, Takara Bio, Cat. No. 630414)) is used to select for the prey library plasmids. SD/-Leu dropout (DO) is so called because the medium includes every essential amino acid except for leucine, which is omitted from the formulation (or “Dropped Out”). Cells harboring Matchmaker prey plasmids are able to grow because pGADT7-Rec encodes the *LEU2* leucine biosynthesis gene, which is otherwise absent from the cell.

## IV. List of Abbreviations

Abbreviation	Definition
AD fusion library (or AD library)	A cDNA library (such as a Mate & Plate Library) constructed in an activation domain (AD) vector such that the proteins encoded by the inserts are fused to the 3' end of the Gal4 AD
AD/library plasmid	Plasmid encoding a fusion of the Gal4 activation domain and a library cDNA
AD/library protein	A protein fusion comprised of the Gal4 activation domain and a polypeptide encoded by a library cDNA
AD vector	Plasmid encoding the yeast Gal4 activation domain

### Yeast Phenotypes

Ade <sup>-</sup> , or His <sup>-</sup> , or Leu <sup>-</sup> , or Trp <sup>-</sup>	Requires adenine (Ade), or histidine (His) or leucine (Leu), or tryptophan (Trp) in the medium to grow, i.e., is auxotrophic for one (or more) of these specific nutrients
LacZ <sup>+</sup>	Expresses the <i>LacZ</i> reporter gene, i.e., is positive for β-galactosidase (β-gal) activity.
Mel1 <sup>+</sup>	Expresses the <i>MEL1</i> reporter gene, i.e., is positive for α-galactosidase (α-gal) activity

### Miscellaneous

SD	Minimal, synthetically defined medium for yeast; is comprised of a nitrogen base, a carbon source (glucose unless stated otherwise), and a DO supplement
DO	Dropout (supplement or solution); a mixture of specific amino acids and nucleosides used to supplement SD base to make SD medium; DO solutions are missing one or more of the nutrients required by untransformed yeast to grow on SD medium. For example, SD/-Leu lacks the amino acid leucine.
YPDA	YPD medium supplemented with adenine hemisulfate (1X concentration = 120 μg/ml)

## V. Control Experiments

**IMPORTANT: Please read the entire protocol before starting.** Use this procedure to perform a control homologous recombination-mediated cloning in yeast BEFORE constructing and screening a two-hybrid library.

### A. General Considerations

The following protocol is designed to familiarize yourself with the procedures and expected results of homologous recombination mediated cloning in Yeast. Additionally, it provides you with a positive control strain that you will need for subsequent yeast two-hybrid screening.

### B. Protocol: Control for Homologous Recombination-Mediated Cloning Protocol

#### 1. Materials:

- Yeastmaker Yeast Transformation System 2
- Y187 Yeast Strain
- pGADT7-Rec Cloning Vector (SmaI-linearized; 500 ng/μl)
- SV40 Large-T PCR fragment (25 ng/μl)
- SD/-Leu with Agar

2. Prepare:
  - Competent Y187 yeast cells, using the Yeastmaker Yeast Transformation System 2 according to the protocol in the accompanying user manual.
  - SD/-Leu agar plates (5 x 100 mm plates; see Appendix C).
3. Combine the following DNA samples in two sterile microcentrifuge tubes and use to cotransform competent Y187 cells:
  - 1 µl pGADT7-Rec (SmaI Linearized)
  - 1 µl pGADT7-Rec (SmaI Linearized) + 3 µl SV40 large T PCR Fragment (25 ng/µl)
4. Plate 100 µl of a 1/10 dilution and 100 µl of a 1/100 dilution on SD/-Leu agar plates.
5. Incubate at 30°C for 3–5 days

**Expected results:** You should see >10X more colonies when cotransforming pGADT7-Rec together with the SV40 Large-T PCR fragment compared to the vector alone. This implies that >90% of the colonies contain a pGADT7-T control plasmid created by homologous recombination.



## VI. Generating the cDNA for Your Library


**IMPORTANT:** Please read the entire protocol before starting. Detailed instructions are provided for first-strand cDNA synthesis (Section A), cDNA amplification using long distance PCR (LD-PCR) (Section B), and column purification of double-stranded (ds) cDNA using a CHROMA SPIN+TE-400 column (Section C).

Use the following protocol for generating cDNA using our simple and highly efficient SMART technology (Figure 2). For a detailed description of SMART technology, refer to Appendix F.


**1 Oligo dT binds to RNA template**



**2 Reverse Transcription**

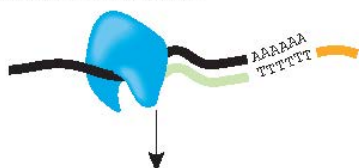

**3 Terminal Transferase Activity adds dCTP; SMART oligo anneals to CCC overhangs**

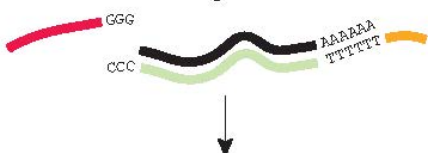

**4 Template Switching incorporates SMART oligo, resulting in known sequences at both ends of all cDNAs**



**Figure 2. Mate & Plate Library construction overview.** SMART cDNA synthesis generates cDNA ends with homology to pGADT7-Rec.

### A. Protocol: First Strand cDNA Synthesis

It is strongly recommended that you perform a positive control cDNA synthesis reaction with Mouse Liver Control Poly A<sup>+</sup> RNA. This control verifies that all components are working properly and lets you compare to the yield and size range of the ds cDNA synthesized from your experimental RNA sample.

**IMPORTANT:** Do not increase the size (volume) of any of the reactions. All components have been optimized for the volumes specified.

The procedure consists of three steps:

- First-strand cDNA synthesis
- Amplification of cDNA by long distance PCR (LD-PCR)
- Column purification of ds cDNA with a CHROMA SPIN+TE-400 column

In the protocol that follows, you have the option of priming first-strand cDNA synthesis with an oligo-dT (CDSIII) or random primer (CDSIII/6). The reaction conditions vary slightly (at Step 6), depending on the primer used.

1. Prepare: high-quality Poly A<sup>+</sup> or total RNA.

**NOTE:** We recommend NucleoSpin RNA kits (Takara Bio, Cat. Nos. 740955.50, 740955.250, or 740955.240C) for purification of total RNA from a variety of sources.

2. Combine and mix the following reagents in a sterile microcentrifuge tube:

1–2 µl	RNA sample (0.025–1.0 µg poly A <sup>+</sup> or 0.10–2.0 µg total RNA)
1.0 µl	CDS III* or CDSIII/6 Primer†
1–2 µl	Deionized H <sub>2</sub> O (to bring volume up to 4.0 µl)
<hr/>	
4.0 µl	Total volume

\*CDSIII = Oligo-dT Primer

†CDSIII/6 = Random Primer

**NOTE:** For the control reaction, use 1 µl (1 µg) of the Mouse Liver Control Poly A<sup>+</sup> RNA .

3. Incubate at 72°C for 2 min
4. Cool on ice for 2 min only, then spin down briefly at 14,000g for 10 sec.
5. Prepare a reaction mix containing the following components, then add this mix to the denatured, primed RNA sample from Step 4, and mix by tapping:

2.0 µl	5X First-Strand Buffer
1.0 µl	DTT (100 mM)
1.0 µl	dNTP Mix (10 mM)
1.0 µl	SMART MMLV Reverse Transcriptase
<hr/>	
9.0 µl	Total volume

6. **ONLY if using Random Primer (CDSIII/6):** Incubate at 25°C for 10 min.

**If using Oligo-dT (CDSIII),** skip this step and continue to Step 7.

7. Incubate at 42° for 10 min.

**NOTE:** Carry out the incubation in a hot lid thermal cycler. If you are using a water bath or non-hot-lid cycler, add a drop of mineral oil to prevent loss in volume due to evaporation.

8. Add 1 µl SMART III-modified oligo, mix and incubate at 42°C for 1 hr.
9. Place the tube at 75°C for 10 min to terminate first-strand synthesis.
10. Cool to room temperature, add 1 µl RNase H (2 units).
11. Incubate at 37° for 20 min.
12. Proceed to LD-PCR amplification (Section VI.B).

**SAFE STOPPING POINT:** Any first-strand synthesis reaction that is not used immediately can be stored at –20°C for up to 3 months.

## **B. Protocol: Amplify cDNA Using Long Distance PCR (LD-PCR)**

Table 2 shows the optimal number of thermal cycles to use based on the amount of RNA used in the first-strand synthesis. Use the minimum number of cycles required to generate 3–6 µg of dsDNA. The optimal cycling parameters were determined using the Mouse Liver Control Poly A<sup>+</sup> RNA; these parameters may vary with different templates and thermal cyclers. This procedure has been tested and optimized only for Advantage 2 Polymerase Mix (Cat. No. 639201).

**Table 2. Relationship between amount of RNA and optimal number of thermal cycles.**

<b>Total RNA (µg)</b>	<b>Poly A<sup>+</sup> RNA (µg)</b>	<b>Number of Cycles</b>
1.0–2.0	0.5–1.0	15–20
0.5–1.0	0.25–0.5	20–22
0.25–0.5	0.125–0.25	22–24
0.05–0.25	0.025–0.125	24–26

- Prepare:
  - First-strand cDNA (results of Section VI.A)
  - Preheat thermal cycler
- Set up TWO 100 µl PCR reactions for each experimental sample and one reaction for the Mouse Liver Control Poly A<sup>+</sup> RNA:

2 µl	First-Strand cDNA (from Section VI.A)
70 µl	Deionized H <sub>2</sub> O
10 µl	10X Advantage 2 PCR Buffer
2 µl	50X dNTP Mix
2 µl	5' PCR Primer
2 µl	3' PCR Primer
10 µl	10X Melting Solution
2 µl	50X Advantage 2 Polymerase Mix
<hr/>	
100 µl	Total volume

- Begin thermal cycling using the following parameters:

95°C    30 sec

N\* cycles:

95°C	10 sec	
68°C	6 min†	

68°C    5 min

\*Refer to Table 2 to estimate the number of cycles, N.


†Program the cycler to increase the extension time by 5 sec with each successive cycle. For example, in the second cycle, the extension should last 6 min and 5 sec; in the third, 6 min and 10 sec, and so on.

4. Analyze a 7  $\mu$ l aliquot of the PCR product from each sample alongside 0.25  $\mu$ g of a 1 kb ladder DNA molecular weight marker on a 1.2% agarose/EtBr gel.

Typical results obtained with Mouse Liver Control Poly A<sup>+</sup> RNA after column chromatography are shown in Figure 3.

**NOTE:** If your PCR product does not appear as expected, refer to the Troubleshooting Guide (Appendix A).

5. Proceed with Column Chromatography (Section VI.C) or store ds cDNA at  $-20^{\circ}\text{C}$  until use.

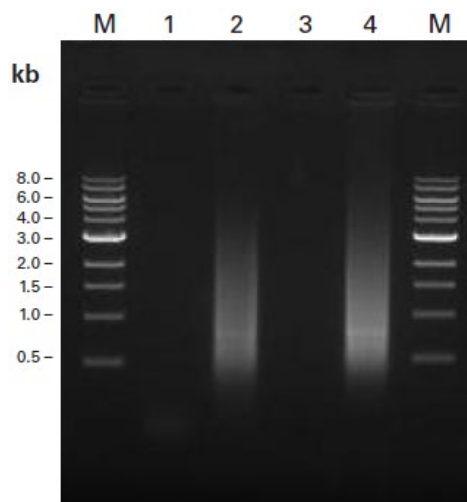

**Figure 3. High-quality cDNA generated using SMART cDNA synthesis.** Oligo dT-primed cDNAs were generated using the Make Your Own “Mate & Plate” Library System. cDNA synthesis was carried out with or without 1  $\mu$ g of Mouse Liver Control Poly A<sup>+</sup> RNA (positive and negative controls, respectively). LD PCR was performed using the Advantage 2 Polymerase Mix (with duplicate samples) and one set of products was purified (size-selected) using CHROMA SPIN+TE-400 columns. Analysis of each respective sample on a 1% agarose gel revealed that the resulting cDNAs ranged from 300 bp to 6 kb. **Lanes M:** 1 kb ladder molecular weight standard. **Lane 1:** unpurified negative control. **Lane 2:** unpurified positive control. **Lane 3:** purified negative control. **Lane 4:** purified positive control. Lane 4 shows reduced abundance of cDNA below 400 bp compared to Lane 2, after size selection with CHROMA SPIN+TE-400 columns.

### C. Protocol: Purify ds cDNA with CHROMA SPIN+TE-400 Columns


In the following protocol, a CHROMA SPIN+TE-400 Column is used to select for DNA molecules  $>200$  bp.

CHROMA SPIN columns are packed with resin that fractionates molecules based on size. Molecules larger than the pore size are excluded from the resin. These molecules quickly move through the gel bed when the column is centrifuged, while molecules smaller than the pore size are retained and trapped inside the gel matrix. For more information about CHROMA SPIN columns, please refer to the [CHROMA SPIN Columns User Manual](#), available at [takarabio.com](http://takarabio.com).

1. Prepare:
  - ds cDNA by LD-PCR (results of Section VI.B; if you do not have at least 3  $\mu$ g, repeat amplification step with more cycles)
  - sodium acetate (3 M)
  - ice cold ethanol (95–100%)

2. Prepare one CHROMA SPIN+TE-400 column for each 93  $\mu$ l cDNA sample (see Figure 4).
  - Invert each column several times to resuspend the gel matrix completely, until it is homogeneous.
  - Remove the top cap.
  - Snap off the break away from the bottom of the column.
  - Place the column in a 2 ml collection tube (supplied).

**NOTE:** You will use 2 columns for each library to be constructed.

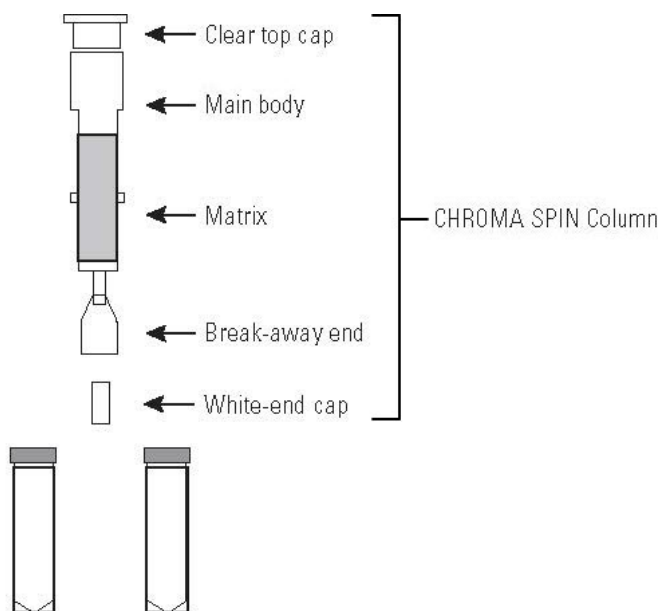

2-ml Collection Tubes

**Figure 4. CHROMA SPIN column and collection tubes.** Note that a conventional, tapered 1.5-ml microcentrifuge tube can be substituted for the 2-ml collection tube when collecting the purified DNA. This will allow the sample to be confined to a narrower area for easier handling. However, a 2-ml collection tube must be used when purging the equilibration buffer.

3. Centrifuge at 700g for 5 min to purge the equilibration buffer, then discard collection tube and buffer. The matrix will appear semi-dry.

**NOTE:** We recommend swing bucket or horizontal rotors. Fixed angle rotors can be used but there is a risk that the sample will pass down the inner side of the columns instead of through the gel matrix, resulting in inconsistent purification.

4. Replace spin column in second collection tube and apply your 93  $\mu$ l sample to the CENTER of the flat surface of the gel matrix.

**NOTE:** Do not allow sample to flow along inner wall of the column.

5. Centrifuge at 700g for 5 min, your purified sample is now in the collection tube.
6. Combine your two purified samples into a single microcentrifuge tube and ethanol-precipitate the cDNA:
  - Add 1/10th vol 3 M sodium acetate.
  - Add 2.5 vol of ice-cold ethanol (95–100%).
  - Place in  $-20^{\circ}\text{C}$  freezer for 1 hr.
  - Centrifuge at 14,000 rpm for 20 min at room temperature.
  - Discard the supernatant; do not disturb the pellet.
  - Centrifuge briefly at 14,000 rpm and remove remaining supernatant.
  - Air dry the pellet for 10 min.

7. Resuspend the cDNA in 20 µl deionized water.

The cDNA is now ready for library construction by in vivo recombination in yeast (Section VII).

**NOTE:** At this point you should have 2–5 µg of ds cDNA.

## VII. Construction of a Two-Hybrid Library

**IMPORTANT:** Please read the entire protocol before starting. Detailed instructions are provided for constructing a “Mate & Plate” yeast two-hybrid library (Section B).

### A. General Considerations for Two-Hybrid Library Construction & Screening

There are two ways to use this kit to make and screen a Matchmaker Two-Hybrid Library.

- The advantage of the recommended protocol is that it provides frozen stocks of the library to screen multiple baits.
- The alternative protocol will not result in a frozen library stock but is a quicker procedure. You can follow this protocol if you do not plan to screen a library multiple times.

#### 1. Recommended Protocol

Make a Mate & Plate library in Y187 Yeast Strain (see Section B, below) and screen against the bait by yeast mating. You will note that the mating procedure is very easy—simply mix a concentrated bait culture with 1 ml of your Mate & Plate library constructed using the protocol and incubate overnight before plating on SD-Leu/-Trp/X- $\alpha$ -Gal/AbA (Aureobasidin A) selective media.

See the [Matchmaker Gold Yeast Two-Hybrid System User Manual](#) for instructions on screening a library using yeast mating.

#### 2. Alternative Protocol

Simultaneously make and screen a library by library-scale cotransformation of bait and prey. Prepare ds cDNA according to the protocol in Section VI and cotransform together with your bait directly into Y2HGold competent cells. Visit [www.takarabio.com](http://www.takarabio.com) for a support protocol.

### B. Protocol: Creating a Mate & Plate Library

#### 1. Prepare:

- Competent Y187 yeast cells using the Yeastmaker Yeast Transformation System 2, according to the protocol in the accompanying user manual.
- 20 µl ds cDNA (results of Section VI)— The amount of ds cDNA should be 2–5 µg.
- The following SD Agar Plates (Appendix C, Section B)
  - SD/-Leu with Agar (100 x 150 mm plates)
  - SD/-Leu with Agar (5–10 x 100 mm plates)
- YPDA/25% Glycerol (Freezing Medium; see Appendix C, Section C)
- Sterile glass beads (5 mm)

**NOTE:** Kanamycin sulfate can be added to all yeast media at a final concentration of 50 µg/ml to stop bacterial contamination (see Appendix C, Section D).

- Follow the library-scale transformation protocol in the [Yeastmaker Yeast Transformation System 2 User Manual](#) to cotransform the following into competent Y187 yeast cells.
  - 20 µl ds cDNA (from Section VI)— The amount of ds cDNA should be 2–5 µg
  - 6 µl pGADT7-Rec (0.5 µg/µl)
- As specified in the transformation protocol ([Yeastmaker Yeast Transformation System 2 User Manual](#), Section VI.B, Step 12), resuspend the transformed cells in 15 ml of 0.9% (w/v) NaCl.
- Spread 100 µl of 1/10 and 1/100 dilutions on SD/-Leu 100 mm agar plates, incubate at 30°C for 3–4 days, and determine the **number of independent clones** in your library (see note below). The number of independent clones is an indication of library complexity.

Refer to Appendix A to improve your transformation efficiency.

### NOTE:

**Number of Independent Clones** = No. of cfu/ml on SD/-Leu x resuspension volume (15 ml)

- This value will determine the transformation efficiency.
  - If you have made a library of 2 million independent clones, you will have 133 colonies on the 1/100 dilution plate.
  - If your library does not have >1 million independent clones, refer to the Appendix A to learn how to improve transformation efficiency and the quality of your cDNA.
  - Only proceed to Step 6 if you have determined that you have >1 million independent clones
- Spread the remainder on 150 mm SD/-Leu Plates
    - 150 µl per plate on ~100 plates
    - Incubate at 30°C for 3–4 days
  - Harvest and Pool Transformants from Step 5.
    - Chill plates at 4°C for 3–4 hr.
    - Add 5 ml of Freezing Medium
    - Use sterile glass beads to detach the colonies from the plate.
    - Combine all liquids in a single sterile flask (you should have a combined total of 0.5 L)
    - Using a hemocytometer, estimate the cell density. If the cell density is  $<2 \times 10^7$  per ml, reduce the volume of the suspension by centrifugation.
    - Aliquot your library into several 1 ml aliquots for short-term use and a few 50 ml aliquots for long-term storage. Store at –80°C.
  - Use a single 1 ml aliquot for library screening using the [Matchmaker Gold Yeast Two-Hybrid System User Manual](#).

**NOTE:** When thawing libraries from storage, always reserve a 10 µl aliquot to recalculate the titer and confirm it remains  $>1 \times 10^7$  cfu/ml. Refer to Appendix B, "Library Titering", for this procedure.



## VIII. References

- For general reviews on yeast two-hybrid systems, see Allen et al., 1995; Bartel et al., 1993a, 1993b; Bartel & Fields, 1997; Fields, 1993; Fields & Sternglanz, 1994; Fritz & Green, 1992; Guarente, 1993; Hopkin, 1996; Luban & Goff, 1995; McNabb & Guarente, 1996; Mendelsohn & Brent, 1994.
- An extensive list of Matchmaker System citations can be obtained from our website ([www.takarabio.com](http://www.takarabio.com)).
- For additional two-hybrid references, see <https://www.caister.com/cimb/v/v1/31.pdf> (Golemis, Serebriiskii, and Law 1999) or use PubMed and search under the key word "two-hybrid" (<https://pubmed.ncbi.nlm.nih.gov/?term=two-hybrid>).

Allen, J. B., Walberg, M. W., Edwards, M. C. & Elledge, S. J. Finding prospective partners in the library: the two-hybrid system and phage display find a match. *Trends Biochem. Sci.* **20**, 511–516 (1995).

Barnes, W. M. PCR amplification of up to 35-kb DNA with high fidelity and high yield from  $\lambda$  bacteriophage templates. *Proc. Natl. Acad. Sci. USA* **91**, 2216–2220 (1994).

Bartel, P., Chien, C.-T., Sternglanz, R. & Fields, S. Elimination of false positives that arise in using the two-hybrid system. *Biotechniques* **14**, 920–924 (1993a).

Bartel, P. L. & Fields, S. (eds). *The Yeast Two-Hybrid System*. (Oxford University Press, Inc., 1997).

Bartel, P. L., Chien, C.-T., Sternglanz, R. & Fields, S. Using the two-hybrid system to detect protein-protein interactions. In *Cellular Interactions in Development: A Practical Approach* (Hartley, D. A.) 153–179 (IRL Press at Oxford University Press, 1993b).

Chenchik, A., Diatchenko, L., Chang, C. & Kuchibhatla, S. Great Lengths cDNA Synthesis Kit for high yields of full-length cDNA. *Clontechiques* **9**, 9–12 (1994).

Chenchik, A., Zhu, Y. Y., Diatchenko, L., Li, R. & Hill, J. Generation and use of high-quality cDNA from small amounts of total RNA by SMART PCR. In *Gene Cloning and Analysis by RT-PCR* (Siebert, P. D. & Larrick, J. W.) 305–319 (Eaton Pub Co, 1998).

Fields, S. The Two-Hybrid System to Detect Protein-Protein Interactions. *Methods* **5**, 116–124 (1993).

Fields, S. & Sternglanz, R. The two-hybrid system: an assay for protein-protein interactions. *Trends Genet.* **10**, 286–292 (1994).

Fritz, C. C. & Green, M. R. Fishing for partners. *Curr. Biol.* **2**, 403–405 (1992).

Golemis, E. A., Serebriiskii, I. & Law, S. F. The Yeast Two-Hybrid System: Criteria for Detecting Physiologically Significant Protein-Protein Interactions. *Curr. Issues Mol. Biol.* **1**, 31–46 (1999).

Golemis, E. A. *et al.* Interaction Trap/Two-Hybrid System to Identify Interacting Proteins. *Curr. Protoc. Mol. Biol.* 20.0.1–20.1.28 (1996).

Guarente, L. Strategies for the identification of interacting proteins. *Proc. Natl. Acad. Sci. USA* **90**, 1639–1641 (1993).

Guthrie, C. & Fink, G. R. (eds). *Guide to Yeast Genetics and Molecular Biology. Methods Enzymol.* **194**, (Academic Press, 1991).

Harper, J. W., Adami, G. R., Wei, N., Keyomarsi, K. & Elledge, S. J. The p21 Cdk-interacting protein Cip1 is a potent inhibitor of G1 cyclin-dependent kinases. *Cell* **75**, 805–816 (1993).

Hopkin, K. Yeast two-hybrid systems: more than bait and fish. *J. NIH Res* **8**, 27–29 (1996).



- James, P., Halladay, J. & Craig, E. A. Genomic Libraries and a Host Strain Designed for Highly Efficient Two-Hybrid Selection in Yeast. *Genetics* **144**, 1425–1436 (1996).
- Kellogg, D. E. *et al.* TaqStart Antibody: ‘hot start’ PCR facilitated by a neutralizing monoclonal antibody directed against Taq DNA polymerase. *Biotechniques* **16**, 1134–1137 (1994).
- Luban, J. & Goff, S. P. The yeast two-hybrid system for studying protein-protein interactions. *Curr. Opin. Biotechnol.* **6**, 59–64 (1995).
- McNabb, D. S. & Guarente, L. Genetic and biochemical probes for protein-protein interactions. *Curr. Opin. Biotechnol.* **7**, 554–559 (1996).
- Mendelsohn, A. R. & Brent, R. Applications of interaction traps/two-hybrid systems to biotechnology research. *Curr. Opin. Biotechnol.* **5**, 482–486 (1994).
- Zhu, Y. Y., Machleder, E. M., Chenchik, A., Li, R. & Siebert, P. D. Reverse Transcriptase Template Switching: A SMART™ Approach for Full-Length cDNA Library Construction. *Biotechniques* **30**, 892–897 (2001).

## Appendix A. Troubleshooting Guide

Table 3. Troubleshooting guide

Problem	Cause	Solution
Low yield of ds cDNA	One or more essential reagents may have been inadvertently omitted from the first-strand or ds cDNA synthesis steps.	Repeat both steps, being careful to check off every item as you add it to the reaction.
	Low yields of ds cDNA may be due to PCR undercycling.	<ul style="list-style-type: none"> <li>• If you suspect undercycling, incubate the PCR-reaction mixture for two more cycles and recheck the product.</li> <li>• If you already used the maximum recommended number of cycles, increase by three more cycles.</li> <li>• If additional cycles do not improve the yield of PCR product, repeat the PCR using a fresh aliquot of the first-strand product.</li> </ul>
	<p>A low yield of ds cDNA may also be due to a low yield of first-strand cDNA. Possible problems with the first-strand reaction include:</p> <ul style="list-style-type: none"> <li>• A mistake in the procedure (such as using a suboptimal incubation temperature or omitting a component), or insufficient RNA in the reaction.</li> <li>• It is also possible that the RNA has been partially degraded (by contaminating RNases) before the first-strand synthesis</li> </ul>	<ul style="list-style-type: none"> <li>• Problems with the first-strand cDNA synthesis can be more easily diagnosed if you perform parallel reactions using the control RNA provided in the kit.</li> <li>• If good results were obtained with the control RNA but not with your experimental RNA, then there may be a problem with your RNA.</li> </ul>

<b>Problem</b>	<b>Cause</b>	<b>Solution</b>
<b>Size distribution of ds cDNA is less than expected</b>	Your RNA starting material may be degraded, very impure, or too dilute.	<ul style="list-style-type: none"> <li>• Check the quality and quantity of your RNA by running a sample on a gel. If the RNA seems too dilute, but otherwise of good quality, restart the experiment using more RNA.</li> <li>• If the RNA seems degraded, restart the experiment using a fresh lot or preparation of RNA.</li> <li>• Also, check the stability of your RNA by incubating a small sample at 37°C for 2 hr. Run it on a gel, parallel to a fresh (unincubated) sample. If the RNA appears to be unstable, it will yield poor results. If this is the case, re-isolate the RNA using a different method (see Related Products).</li> <li>• Problems with your RNA are easily diagnosed if you perform parallel reactions using the control RNA provided in this kit.</li> </ul>
<b>Presence of low molecular weight (&lt;0.1 kb) material in the ds cDNA product</b>	<ul style="list-style-type: none"> <li>• The raw cDNA (e.g., before size fractionation) is expected to contain some low-molecular weight DNA contaminants, including unincorporated primers, SMART oligonucleotides, and very short PCR products.</li> <li>• However, these small fragments are generally removed from the ds cDNA preparation in the size fractionation step using the columns provided.</li> <li>• Note that a preponderance of low-molecular-weight (&lt;0.1 kb) material in the raw PCR product may be indicative of overcycling.</li> </ul>	If you suspect overcycling, then the PCR step must be repeated with a fresh sample of first-strand cDNA, using 2–3 fewer cycles.

<b>Problem</b>	<b>Cause</b>	<b>Solution</b>
<b>Presence of low molecular weight (&lt;0.1 kb) in the size-fractionated ds cDNA</b>	CHROMA SPIN columns are designed to remove low-molecular weight cDNA fragments, small DNA contaminants, and unincorporated nucleotides from the cDNA.	The column should be stored and used at room temperature. If it is chilled at 4°C and then warmed to room temperature for use, bubbles may form, which interfere with the proper functioning of the column.
	The resolving function of the column will be diminished if the gel matrix becomes dry. In drying, the matrix body may shrink away from the innerwall of the column casing. The ds cDNA mixture can then flow down the sides of the column, allowing small contaminants to elute with the cDNA.	Extreme, uneven deposition of the ds cDNA mixture on the surface of the column can cause inefficient separation of ds cDNA from low-molecular-weight contaminants.
<b>Low transformation efficiency</b>	DNA is not sufficiently pure	Check the purity of the DNA and, if necessary, repurify by ethanol precipitation.
	The fusion protein encoded by the DNA-BD/bait plasmid (Two-Hybrid System) may be toxic.	Try using a vector that expresses lower levels of the fusion, such as pGBT9.
	Improper media preparation	Remake media and test with control transformations. Check the efficiency using the pGBT9 Control Plasmid. Plate on SD/-Trp. The transformation should yield $\geq 1 \times 10^5$ colonies/ $\mu$ g DNA.
	Slow-growing Y187 yeast cells	Streak to single colonies, inoculate 4 large colonies into YPDA broth, and select the fastest growing to make competent cells.
<b>Yeast two-hybrid screening problems</b>		See <a href="http://www.takarabio.com">Matchmaker Gold Yeast Two-Hybrid System User Manual</a> at <a href="http://www.takarabio.com">www.takarabio.com</a>

## Appendix B. Library Titering

### A. General Considerations

- Diluted libraries are always less stable than undiluted libraries. Therefore, once the library dilutions are made, plate them within the next hour, before misleading reductions in titer can occur.
- Use proper sterile technique when aliquoting and handling libraries.
- Design and use appropriate controls to test for cross-contamination.
- Always use the recommended concentration of antibiotic (i.e., kanamycin) in the medium to suppress growth of contaminating bacteria.

## B. Library Titering

Materials:

- YPDA Broth (Appendix C)
- SD/-Leu (100 mm plates) (Appendix C)
- Sterile glass spreading rod, bent Pasteur pipette, or 5 mm glass beads for spreading culture on plates

1. Transfer the 10 µl library aliquot (reserved from Section VII.B) to 1 ml of YPDA Broth in a 1.5 ml microcentrifuge tube. Mix by gentle vortexing. This is **Dilution A** (dilution factor = 10<sup>-2</sup>).

**NOTE:** Matchmaker Gold Libraries are extremely viscous, so pay special attention during aliquoting.

2. Remove 10 µl from **Dilution A** and add it to 1 ml of YPDA Broth in a 1.5 ml microcentrifuge tube. Mix by gentle vortexing. This is **Dilution B** (dilution factor = 10<sup>-4</sup>).
3. Add 10 µl from **Dilution A** to 50 µl of YPDA Broth in a 1.5-ml microcentrifuge tube. Mix by gentle vortexing. Spread the entire mixture onto an SD/-Leu plate.
4. Remove 50 µl aliquots from **Dilution B** and spread onto separate SD/-Leu plates as above.
5. Invert the plates and incubate at 30°C for 3–5 days.

**NOTE:** Colony size will vary, depending on the insert.

6. Count the number of colonies on plates having 30–300 colonies.
7. Calculate the titer (cfu/ml) as follows:

$$\frac{\text{Number of colonies}}{\text{plating volume (ml) x dilution factor}} = \text{cfu/ml}$$

### NOTES:

- Due to slight variability in pipettes and pipetting techniques, a 2–5-fold range in titer calculations is not unusual.
- Example calculation:
  - No. of colonies on plate = 100
  - Plating volume = 0.05 ml
  - Dilution factor = 10<sup>-4</sup>

$$\frac{100}{0.05 \text{ ml} \times 10^{-4}} = 2 \times 10^7 \text{ cfu/ml}$$

## Appendix C. Yeast Growth Media & Supplements

### A. Ready-to-go Media Pouches Available from Takara Bio

We offer media sets with a complete assortment of mixes in convenient, “ready-mixed” foil pouches, for use with the Matchmaker Gold Yeast Two-Hybrid System.

- See Table 4 for a list of the components of the Yeast Media Set 2 (Cat. No. 630494) and the Yeast Media Set 2 Plus (Cat. No. 630495).

**NOTE:** The Yeast Media Set 2 Plus also contains the additional media supplements Aureobasidin A and X-α-Gal, which are required for yeast two-hybrid screening using the Matchmaker Gold Yeast Two-Hybrid System, but not with this user manual.

## Make Your Own “Mate & Plate™” Library System User Manual

- See Table 5 for information for purchasing each of the media mixes separately, in packs of 10 pouches.
- See Table 6 for information on preparing all additional required media supplements and purchasing Aureobasidin A and X- $\alpha$ -Gal separately.

**Table 4. Components of Yeast Media Set 2 & Yeast Media Set 2 Plus.**

Media Pouch	Number of Pouches	Volume of Media (per Pouch)
YPDA Broth	2	0.5 L
YPDA with Agar	1	0.5 L
SD/-Leu Broth	1	0.5 L
SD/-Leu with Agar	1	0.5 L
SD/-Trp Broth	1	0.5 L
SD/-Trp with Agar	1	0.5 L
SD/-Leu/-Trp with Agar	10	0.5 L
SD/-Ade/-His/-Leu/-Trp with Agar	1	0.5 L

### Additional Components in Yeast Media Set 2 Plus

X- $\alpha$ -Gal	250 mg	—
Aureobasidin A	1 mg	—

**Table 5. Individual Yeast Media Pouches for Matchmaker Gold Protocols.**

Yeast Media Pouches	Cat. No.	Number of Pouches	Volume of Media (per Pouch)
<b>Rich Media (for Routine Culturing of Untransformed Yeast)</b>			
YPDA Broth	630306	10	0.5 L
YPDA with Agar	630307	10	0.5 L
<b>Minimal Media Single Dropouts (SDO)</b>			
SD/-Trp Broth	630308	10	0.5 L
SD/-Trp with Agar	630309	10	0.5 L
SD/-Leu Broth	630310	10	0.5 L
SD/-Leu with Agar	630311	10	0.5 L
<b>Minimal Media Double Dropouts (DDO)</b>			
SD/-Leu/-Trp Broth	630316	10	0.5 L
SD/-Leu/-Trp with Agar	630317	10	0.5 L
<b>Minimal Media Triple Dropouts (TDO)</b>			
SD/-His/-Leu/-Trp Broth	630318	10	0.5 L
SD/-His/-Leu/-Trp with Agar	630319	10	0.5 L
<b>Minimal Media Quadruple Dropouts (QDO)</b>			
SD/-Ade/-His/-Leu/-Trp Broth	630322	10	0.5 L
SD/-Ade/-His/-Leu/-Trp with Agar	630323	10	0.5 L

**Table 6. Additional Media and Media Supplements Required for a Two-Hybrid Screen.**

Freezing Medium		Preparation Instructions
YPDA Medium & 25% glycerol		see Appendix C, Section C

Supplement Name	Cat. No.	Stock Solution Concentration
Aureobasidin A	630466	500 µg/ml
X-α-Gal (250 mg)	630463	20 mg/ml in dimethyl formamide
Kanamycin sulfate	—	50 mg/ml stock solution
Dimethyl formamide	—	—

## B. General Media Preparation Instructions

- Prepare media by dissolving pouch contents in 500 ml ddH<sub>2</sub>O, autoclave for 15 min at 121°C, and allow to cool before use (or filter-sterilize broth media). Do not over-autoclave.
- This media does not usually require pH adjustment, but if your source water is particularly acidic, you may need to adjust the pH of the media to 5.8.
- For additional information on preparing media from the pouches, please see the [Yeast Media Protocol-at-a-Glance](#) at [www.takarabio.com](http://www.takarabio.com)

## C. Freezing Medium

Mix 100 ml YPDA (sterile) and 50 ml 75% glycerol (sterile).

## D. Kanamycin Supplement to Yeast Media

Kanamycin sulfate can be added to all yeast media at a final concentration of 50 µg/ml to stop bacterial contamination. Please note that kanamycin does not stop contaminating fungal growth, so proper sterile technique must still be used. Also note that this antibiotic does not select for any plasmids in yeast.

## Appendix D. Host Strain Information

The phenotype and complete genotype of Y187 Yeast Strain (the library strain) are shown in Tables 7 and 8. For additional information on the growth and maintenance of yeast, see the supporting Matchmaker protocols at [www.takarabio.com](http://www.takarabio.com). We also recommend the Guide to Yeast Genetics and Molecular Biology (Guthrie & Fink, 1991).

**Table 7. Yeast Host Strain Genotype.**

Strain	Genotype*	Reporters	Transformation Markers	Reference
Y187 <sup>†</sup>	<i>MATα, ura3-52, his3-200, ade2-101, trp1-901, leu2-3, 112, gal4Δ, gal80Δ, met-</i> , <i>URA3 :: GAL1<sub>UAS</sub>-Gal1<sub>TATA</sub>-LacZMEL1</i>	<i>MEL1, LacZ</i>	<i>trp1, leu2</i>	Harper et al., 1993


\*The *GAL1*, *GAL2*, and *MEL1* upstream activating sequences (UASs) are recognized and bound by the Gal4 BD. The *trp1*, *his3*, *gal4*, and *gal80* mutations are all deletions; *leu2-3, 112* is a double mutation.

<sup>†</sup>The *LacZ* reporter construct was integrated into the yeast genome by homologous recombination at the *ura3-52* mutation. Recombinants were selected on SD/-Ura.


**Table 8. Phenotype Testing on Various SD Media.**

Strain	SD/-Ade	SD/-His	SD/-Leu	SD/-Trp	SD/-Ura
Y187	–	–	–	–	+
Y187 (pGADT7-T)	–	–	+	–	+

Appendix E. Plasmid Information

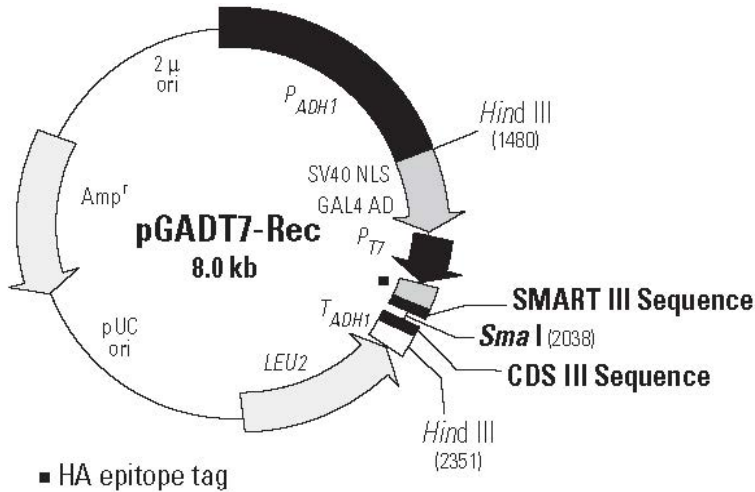

SMART III terminus

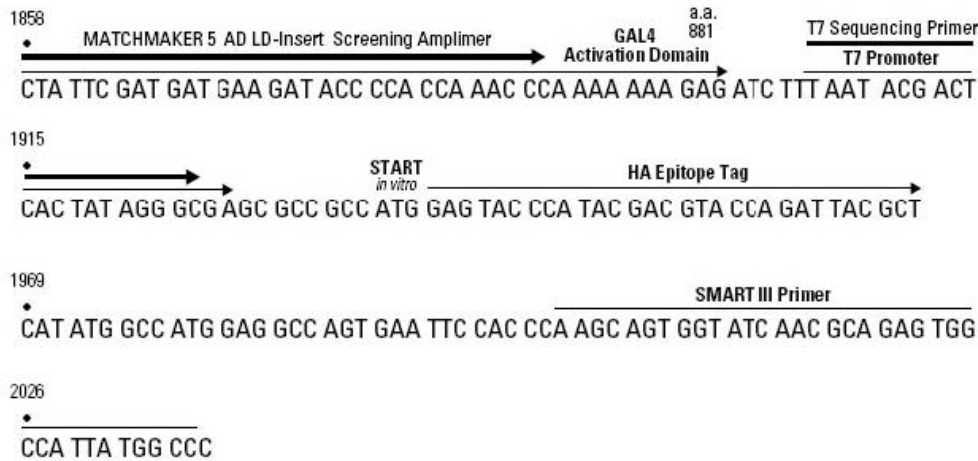

CDS III terminus

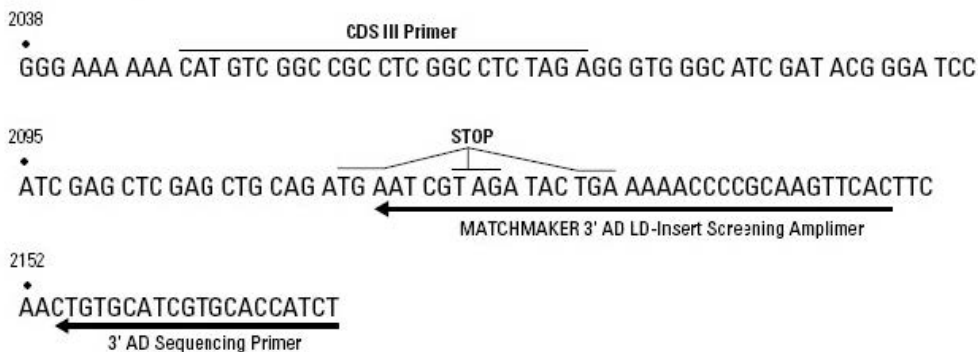


Figure 5. Restriction Map and Cloning Sites of pGADT7-Rec. A unique restriction site (SmaI) is shown in bold. The Make Your Own "Mate & Plate™" Library System (Cat. No. 630490) contains the SmaI-linearized form of this vector, the form used for recombination-mediated cloning in yeast.



## Appendix F. SMART Technology Overview

### A. SMART Technology

Messenger RNA transcripts are efficiently copied into ds cDNA using SMART (Switching Mechanism at 5' end of RNA Template) technology (Zhu, Y.Y., et al., 2001). This cDNA synthesis and amplification system is particularly well suited for two-hybrid library construction because it consistently delivers high yields of cDNA while maintaining sequence representation. By maintaining the complexity of the original tissue, the SMART procedure provides you with the best opportunity to detect rare and potentially novel interactions during yeast two-hybrid screening.

### B. Mechanism of cDNA Synthesis

In the first-strand cDNA synthesis step, MMLV (Moloney Murine Leukemia Virus) Reverse Transcriptase (RT) is used to transcribe RNA into DNA. To prime RNA for cDNA synthesis, you may use either a modified oligo(dT) primer (our CDS III Primer) or a random primer (our CDS III/6 Primer).

The composition of the resulting cDNA library may differ depending on which primer you choose. If you use the CDS III Primer, which hybridizes to the 3'-end of poly A<sup>+</sup> RNA, sequences close to the 5'-end of the transcript may be slightly under-represented. If instead you use the CDS III/6 Primer, a random primer that can hybridize to many different sequences on the RNA template, your library should contain a variety of 5'- and 3'-end sequences, which are represented in near equal proportions.

When MMLV RT encounters a 5'-terminus on the template, the enzyme's terminal transferase activity adds a few additional nucleotides, primarily deoxycytidine, to the 3' end of the cDNA. The SMART III Oligonucleotide, which has an oligo(G) sequence at its 3' end, base-pairs with the deoxycytidine stretch, creating an extended template (Figure 2, above). RT then switches templates and continues replicating to the end of the oligonucleotide. In the majority of syntheses, the resulting ss cDNA contains the complete 5' end of the mRNA as well as the sequence complementary to the SMART III Oligo, which then serves as a universal priming site (SMART anchor) in the subsequent amplification by long-distance PCR (LD PCR; Chenchik *et al.*, 1998). Only those ss cDNAs having a SMART anchor sequence at the 5' end can serve as a template and be exponentially amplified by long-distance PCR (LD PCR).

In the second step, ss cDNA is amplified by LD PCR to produce a ds cDNA library. We recommend using the Advantage 2 PCR Polymerase Mix (Cat. No. 639201) to generate and amplify ds cDNA. The Advantage 2 Polymerase Mix consists of Titanium® Taq DNA Polymerase (a nuclease-deficient N-terminal deletion of Taq DNA polymerase), TaqStart® Antibody to provide automatic hot-start PCR (Kellogg et al., 1994), and a minor amount of a proofreading polymerase. This polymerase system lets you amplify cDNA (as large as 20 kb) with a fidelity rate significantly higher than that of conventional PCR (Barnes, 1994).



Contact Us	
Customer Service/Ordering	Technical Support
tel: 800.662.2566 (toll-free)	tel: 800.662.2566 (toll-free)
fax: 800.424.1350 (toll-free)	fax: 800.424.1350 (toll-free)
web: <a href="http://takarabio.com/service">takarabio.com/service</a>	web: <a href="http://takarabio.com/support">takarabio.com/support</a>
e-mail: <a href="mailto:ordersUS@takarabio.com">ordersUS@takarabio.com</a>	e-mail: <a href="mailto:technical_support@takarabio.com">technical_support@takarabio.com</a>

## Notice to Purchaser

Our products are to be used for **Research Use Only**. They may not be used for any other purpose, including, but not limited to, use in humans, therapeutic or diagnostic use, or commercial use of any kind. Our products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products or to provide a service to third parties without our prior written approval.

Your use of this product is also subject to compliance with any applicable licensing requirements described on the product’s web page at [takarabio.com](http://takarabio.com). It is your responsibility to review, understand and adhere to any restrictions imposed by such statements

© 2021 Takara Bio Inc. All Rights Reserved.

All trademarks are the property of Takara Bio Inc. or its affiliate(s) in the U.S. and/or other countries or their respective owners. Certain trademarks may not be registered in all jurisdictions. Additional product, intellectual property, and restricted use information is available at [takarabio.com](http://takarabio.com).

This document has been reviewed and approved by the Quality Department.