
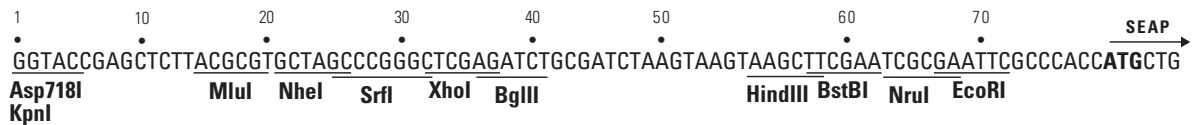


pSEAP2-Basic Vector Information

PT3075-5

GenBank Accession #: U89937

Cat. Nos. 631715 & 631735

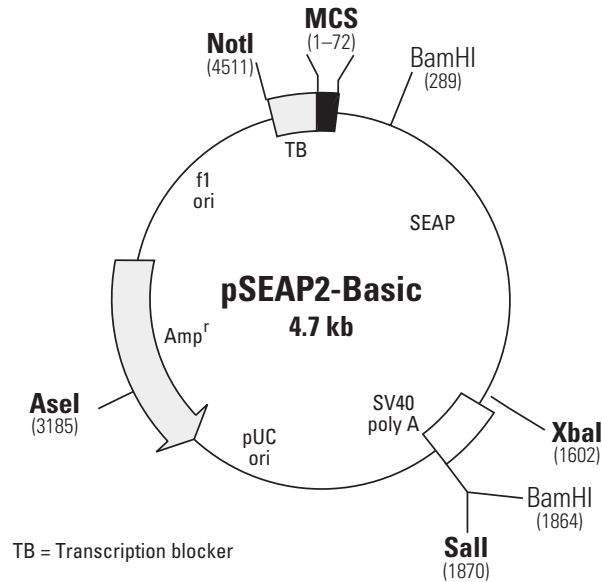

Restriction Map and Multiple Cloning Site (MCS) of pSEAP2-Basic. Unique restriction sites are in bold.

Description:

pSEAP2-Basic allows for expression of the reporter gene secreted alkaline phosphatase (SEAP). This vector lacks eukaryotic promoter and enhancer sequences and has an MCS that allows putative promoter DNA fragments to be inserted upstream of the SEAP gene. Enhancers can be cloned into either the MCS or unique downstream sites. The SEAP coding sequence is followed by the SV40 late polyadenylation signal to ensure proper, efficient processing of the SEAP transcript in eukaryotic cells. A synthetic transcription blocker (TB), composed of adjacent polyadenylation and transcription pause sites, located upstream of the MCS reduces background transcription (1). The vector backbone also contains an f1 origin for single-stranded DNA production, a pUC origin of replication, and an ampicillin resistance gene for propagation and selection in *E. coli*. The SEAP2 Vectors incorporate a number of features that improve the sensitivity of SEAP by increasing the efficiency of SEAP expression or that enhance the utility of the vectors. These include: an improved Kozak consensus translation initiation site (2); the removal of the SV40 small-t intron, which can cause cryptic splicing and reduced expression in some genes and/or cell types (3, 4); switching from the early to late polyadenylation signal of SV40, which typically causes a five-fold increase in mRNA levels (5); an expanded multiple cloning site (MCS); compact plasmid size; and removal of extraneous sequences from the 3' untranslated region of the SEAP mRNA.



Clontech

United States/Canada
800.662.2566

Asia Pacific
+1.650.919.7300

Europe
+33.(0)1.3904.6880

Japan
+81.(0)77.543.6116

Clontech Laboratories, Inc.
A Takara Bio Company
1290 Terra Bella Ave.
Mountain View, CA 94043
Technical Support (US)
E-mail: tech@clontech.com
www.clontech.com

(PR6X2122; published 30 October 2006)

Use:

The integrated set of Great EscAPE™ SEAP2 Vectors has been designed to provide maximal flexibility in studying regulatory sequences from the gene of interest. pSEAP2-Basic allows for expression of SEAP under control of putative promoters and/or enhancers of interest. The secreted SEAP enzyme is assayed directly from the culture medium and permits time-course studies not possible with assays dependent on cell lysates. Furthermore, the cells can be used for further investigations such as RNA or protein studies. The SEAP2 Vectors can be transfected into mammalian cells by any standard method.

Location of features:

- Multiple Cloning Site: 1–72
- Secreted alkaline phosphatase (SEAP) gene:
 - Kozak consensus translation initiation signal: 75–82
 - SEAP coding sequences:
 - Start codon (ATG): 80–82
 - Stop codon: 1637–1639
 - Signal peptide: 80–130
 - Mature protein: 131–1636
 - C-terminal extension to SEAP: 1598–1636
- SV40 Late mRNA polyadenylation signal: 1750–1755
mRNA 3' end: 1769
- pUC plasmid replication origin: 2148–2791
- Ampicillin resistance gene:
 - Promoter: –35 region: 3869–3864; –10 region: 3846–3841
 - Transcription start point: 3834
 - Ribosome binding site: 3811–3807
 - β-lactamase coding sequences:
 - start codon (ATG): 3799–3797; stop codon: 2941–2939
 - β-lactamase signal peptide: 3799–3731
 - β-lactamase mature protein: 3730–2942
- f1 single-strand DNA origin (packages the coding strand of SEAP): 3931–4386
- Transcription blocker (TB): 4517–4670
 - Synthetic polyadenylation site (6): 4517–4565
 - Transcription pause site from human α2 globin gene (7): 4579–4670

Recommended sequencing primers:

- 5' of MCS: 4619–4638 (5'-CTAGCAAATAGGCTGTCCC-3')
- 3' of MCS: 185–165 (5'-CCTCGGCTGCCTCGCGTTCC-3')

Propagation in *E. coli*

- Suitable host strains: DH5α and other general purpose strains. Single-stranded DNA production requires a host containing an F' episome such as JM109.
- Selectable marker: plasmid confers resistance to ampicillin (100 µg/ml) to *E. coli* hosts.
- *E. coli* replication origin: pUC
- Copy number: ~500
- Plasmid incompatibility group: pMB1/Col E1

References:

1. Eggermont, J. & Proudfoot, N. (1993) *EMBO J.* **12**:2539–2548.
2. Kozak, M. (1987) *Nucleic Acids Res.* **15**:8125–8148.
3. Huang, M. T. F. & Gorman, C. M. (1990) *Mol. Cell. Biol.* **10**:1805–1810.
4. Evans, M. J. & Scarpulla, R. C. (1989) *Gene* **84**:135–142.
5. Carswell, S. & Alwine, J. C. (1989) *Mol. Cell. Biol.* **9**:4248–4258.
6. Levitt, N. *et al.* (1989) *Genes Dev.* **3**:1019–1025.
7. Enriquez-Harris, P. *et al.* (1991) *EMBO J.* **10**:1833–1842.

Note: The attached sequence file has been compiled from information in the sequence databases, published literature, and other sources, together with partial sequences obtained by Clontech. This vector has not been completely sequenced.

Notice to Purchaser

This product is intended to be used for research purposes only. It is not to be used for drug or diagnostic purposes, nor is it intended for human use. Clontech products may not be resold, modified for resale, or used to manufacture commercial products without written approval of Clontech Laboratories, Inc.

Clontech, Clontech logo and all other trademarks are the property of Clontech Laboratories, Inc.

Clontech is a Takara Bio Company. ©2006